


Helpful Tips

How to Charge Your Lively Mobile2

When you're not using your Lively Mobile2, always place it in the charging dock so it's ready to use. We recommend charging when you go to bed.

The battery indicator will flash red and a tone will play, when the battery is low. The tone will repeat more frequently as the battery level decreases. If the battery loses all charge, the device will turn off. Once placed on the charging dock, the Lively Mobile2 will power back up.

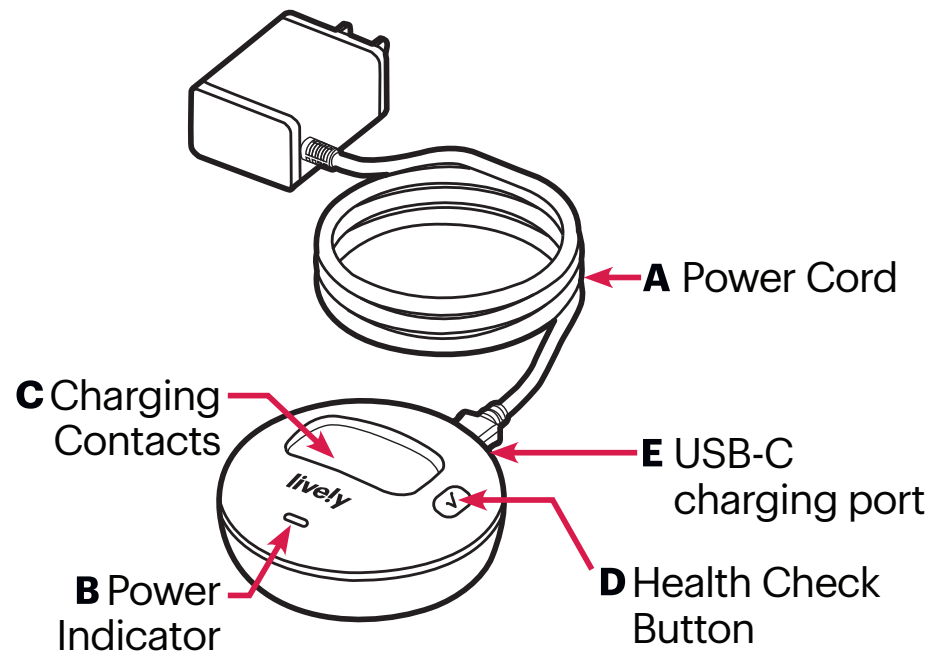

A. Power Cord – Provides power to the charging dock from a wall outlet.

B. Power Indicator – Lights up green when plugged into a wall outlet and blinks green when charging the Lively Mobile2.

C. Charging Contacts – Allows charging of your Lively Mobile2 when docked.

D. Health Check Button – Check that your Lively Mobile2 is operating correctly.

E. USB-C Charging Port - Connects the Power Cord to the Charging Dock